



Taha Kimia Tajhiz Co.




Core Lab

# **Well-Site Gamma Logger, WSGL-100 Datasheet**

**Routine Core Analysis Laboratory (RCAL)/ Core  
Gamma Logging**



## Well-Site Gamma Logger, WSGL-100

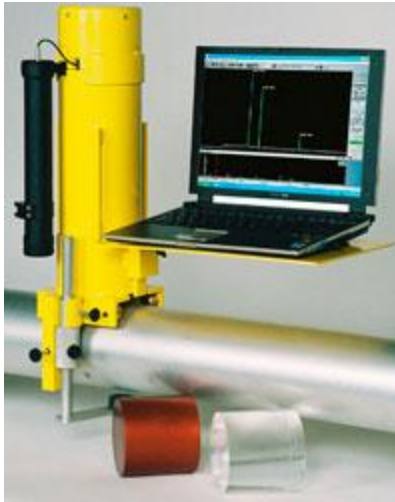

The WSGL-100, Well-site Gamma Logger is a portable gamma logger complete with 3" diameter NaI crystal detector and photomultiplier; lead shielding tunnel; MiNiH smart battery; onboard computer; flight case; battery case; operating software; software for data transfer; 3Kg KCl for system calibration; operating manuals; export packing.

System includes a 5 micro-curie cesium source for spectral calibration and aluminum and polycarbonate standards for bulk density calibration. A Windows based graphical interface easily and quickly guides the operator through the calibration and logging procedures storing data to the hard drive.

